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This manual contains important
informations about this product.
Please pass it to the final user
upon delivery.

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ENGLISH

Physipro inc. is delighted to have you among its clients. We sincerely wish to thank you for the trust you have shown in our Company by choosing one of our products. This user manual has been designed to allow the client to use the Symbio headrest support, a Physipro inc. product, safely in an optimal way. For all adjustments and settings, Physipro inc. strongly recommends you consult a qualified professional.

Physipro inc. provides you with specific verifications you should perform on a regular basis in order to optimize the performance of your Symbio headrest and extend its operational life. Physipro inc. relieves itself from any liability should any bodily injury or property damage result from a lack of care or misuse or from any modification made to the product without prior written consent of the Company.


For Physipro inc., your satisfaction remains a priority.

Table of Contents

1. Components	2
2. Recommendations	3
3. Installation on Rigid Backrest	3
4. Installation on Tension Bar	4
5. Adjusting the Headrest Support	4
6. Lower and Upper Sections Depth Adjustment	4
7. Angle Adjustment	5
8. 45° Axis Adjustment	5
9. Multidirectional Adjustment	5
10. Warranty	6

1. Components

1. Anchoring ring (AT5210)
2. Fixture (PH-02066B)
3. Locking ring (PH-02066C)
4. Lower section (PH-02386)
5. Depth adjustment hinge (PH-02387)
6. Joint for headrest support mechanism (PH-02388)
7. Lateral hinge fitting (PH-02389)
8. Fixing screw M5 (16mm) (ST-51387)
9. Fixing screw M6 pivot (20 mm) (ST-51464)
10. Handle (ST-52034)

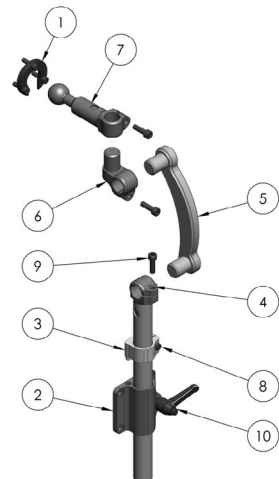2. Recommendations

Do not use this equipment without having read and understood this user manual in its entirety. This manual includes critical information intended to ensure the safety of the user and the individuals who handle the wheelchair.

2.1 Basic verification

At the time of delivery and receipt check all parts for shipping damage before using.


2.2 Safety rules

Several safety measures must be applied to ensure the safety of the user. Note that the following list is not comprehensive. Therefore, the people who use the Symbio headrest have the responsibility to use it cautiously:

- ✓ Do not use the Symbio headrest to push or lift wheelchair.
- ✓ Installation and adjustment of your Symbio headrest should be performed only by an authorized dealer or qualified therapist.
- ✓ Make sure that all screws and adjustments are tightened properly, otherwise serious injury damage may result.

3. Installation on Rigid Backrest


- 3.1 Place the mounting bracket to about 4" (10 cm) from the top of the backrest.
- 3.2 Drill the backrest by using the mounting bracket as a guide.
- 3.3 Install the bracket screws and nuts as shown.

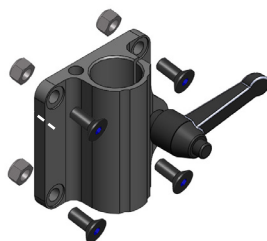ENGLISH


4. Installation on Tension Bar

- 4.1 Align the adapter **A** into the interior of the rod with the bracket **B**.
- 4.2 Install the bracket screws and nuts as shown.


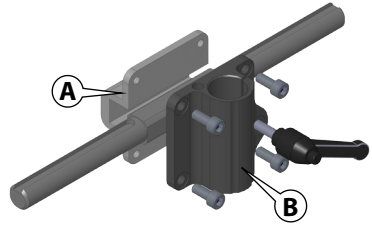
5. Adjusting the Headrest Support

- 5.1 Align headrest tube with the hole in the mounting bracket and slide into mounting bracket on back.
- 5.2 Adjust the headrest hardware to the desired position and secure in place by tightening the seat back mount handle **A**.
- 5.3 To adjust the height, slide the locking ring **B** against the bracket and tighten.

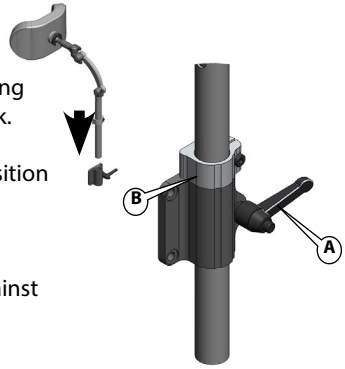
6. Lower and Upper Sections Depth Adjustment

- 6.1 Loosen screws (M6) and rotate the horizontal adjustment joint to let the depth adjustment tube slide freely.
- 6.2 Adjust the headrest to the desired position.

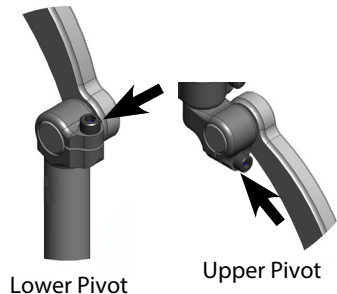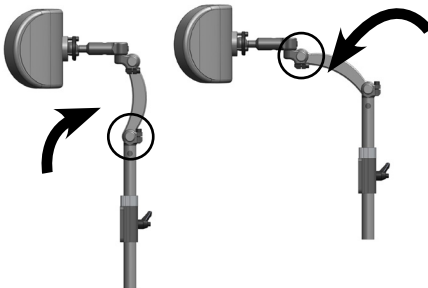
7. Angle Adjustment

- 7.1 Loosen screw (M6) to let the angle adjustment tube slide freely.

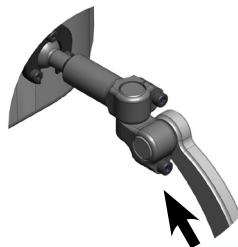- 7.2 Adjust the angle adjustment tube in the desired position and securely tightens screws (M6).

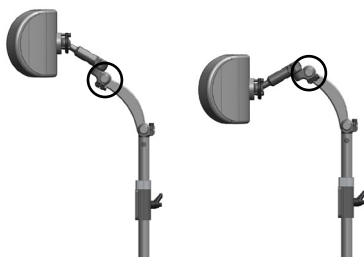8. 45° Axis Adjustment


- 8.1 Loosen screw (M6) to let the axis adjustment tube slide freely.

- 8.2 Adjust the axis in the desired position and securely tightens screws (M6).

9. Multidirectional Adjustment

- 9.1 Loosen screws (M6) to make the adjustment.

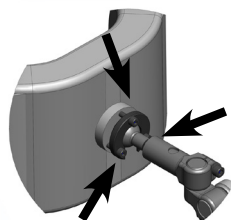ENGLISH

10. Warranty

The Symbio headrest support, a product of Physipro inc., is covered with our twenty-four (24) months manufacturer's warranty. It includes parts and labor on manufacturer's defects and optional parts, and all new parts replaced by your dealer or an authorized establishment from the date of delivery or replacement.

Physipro inc. will repair or replace the defective parts as long as the warranty is valid. You will have to contact Physipro inc. or an authorized dealer to obtain your warranty services. Do not return your backrest without prior consent. If you are not satisfied with the repairs, please send us a written complaint with your dealer's name, address, delivery date and the product's serial number to the address on this manual cover.

Exclusions and limitations

This warranty does not apply to products that were intentionally or accidentally altered, misused, neglected, did not have the proper maintenance or had the serial number removed or erased. Furthermore, this warranty does not apply if repairs or modifications were done without a written consent from Physipro inc., or an authorized dealer. The same restrictions apply if any other kind of damages were done without Physipro inc. involvement. This warranty will not apply for everyday use or failure to respect the rules in this manual.

Physipro inc. will not take any responsibilities for damages incurred during shipping.

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